



REF: 908-L

REF: 408-L

REF: 308-L

- Available with 4 different tips
- Choose between finger or foot controlled instruments
- Less smoke and char
- Cost effective

LiNA Endo-Hooks

Single use endoscopic monopolar hook electrodes for laparoscopic procedures

LiNA Endo-Hooks

The single use LiNA monopolar Endo-Hook provides the right solution for today's demand for laparoscopic procedures. The range of lightweight and slim LiNA Endo-Hooks offers excellent and precise dissection, cutting and coagulation. The ceramic tip insulation leads the energy exactly where it is needed, thus ensuring full safety throughout the procedure. The variety of 4 tip configurations (J-shaped, L-shaped, Needle and Spatula) gives the best function for all surgical applications. Choose between finger and foot-operated instruments.

The standard length of the LiNA 308, 408 and 908 Endo-Hook is 320mm. Furthermore the foot-operated REF 308-L is also available in 400mm length for larger patients, i.e. obesity patients, order REF 308-L-400.

The LiNA Endo-Hook is available in 3 different versions:

308 – Foot-operated without cable

The 308 Endo-Hook is the most cost effective version, where you attach a reusable cable to the 4mm connector in the 308 handle.

LiNA offers a reusable 3m silicone cable REF REC-617 with 4mm generator connector. Adaptor with 8mm generator connector is also available – LiNA REF 572.


408 – Foot-operated with 3m cable with 4mm connector

The range of 408 Endo-Hooks includes a 3m cable with 4mm generator connector. Adaptor with 8mm generator connector is also available – LiNA REF 572.

908 – Finger-operated with 3m cable and US-3 pin connector

The range of 908 Endo-Hooks offers the surgeon a finger-operated Endo-Hook, which eliminates the use of a foot pedal during the procedure. The 908 Endo-Hook is the surgeon's preferred choice because working with the 908 Endo-Hook is like working with a standard finger-operated diathermy pencil.

The 908 Endo-Hook includes a 3m cable and a US-3 pin generator connector. Adaptors for most generators are also available, see below list of REF 571 generator adaptors.

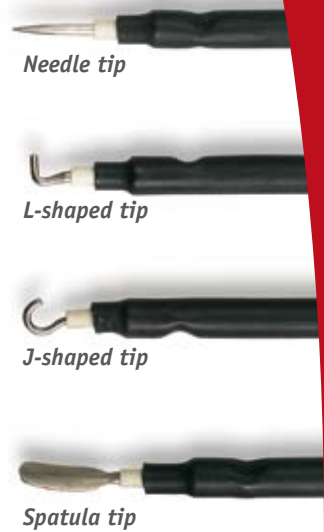

Needle tip

L-shaped tip

J-shaped tip

Spatula tip



REF: 572

REF: 571-

Ordering information:

LiNA Endo-Hooks 0459

REF.	Product	Packing
308-J (J-shaped tip) 308-L (L-shaped tip) 308-N (Needle tip) 308-S (Spatula tip)	LiNA Endo-Hook - 320mm instrument – Foot controlled instrument with 4mm male connector in the handle, excl. cable – Single use. Available with 4 different tips.	10 pcs per box, sterile.
308-L-400 (L-shaped tip)	LiNA Endo-Hook - As above but with extra long 400mm instrument.	10 pcs per box, sterile.
408-J (J-shaped tip) 408-L (L-shaped tip) 408-N (Needle tip) 408-S (Spatula tip)	LiNA Endo-Hook – 320mm instrument – Foot controlled instrument with 3m integrated cable and 4mm generator connector – Single use. Available with 4 different tips.	10 pcs per box, sterile.
908-J (J-shaped tip) 908-L (L-shaped tip) 908-N (Needle tip) 908-S (Spatula tip)	LiNA Endo-Hook – 320mm instrument – Finger controlled instrument with 3m integrated cable and US-3 pin standard generator connector – Single use. Available with 4 different tips.	10 pcs per box, sterile.
572 (for 308 & 408)	LiNA HF generator Adaptor – reusable – 8mm connector for foot controlled endo-hooks. (for Valleylab generators).	1 pc per box.
571-ACC (for 908) 571-AE 571-B 571-E 571-L 571-M 571-N 571-S1 571-S2	LiNA HF generator Adaptor – reusable – For Erbe ACC generator LiNA HF generator Adaptor – reusable – For Aesculab generator LiNA HF generator Adaptor – reusable – For Bovie generator LiNA HF generator Adaptor – reusable – For Old Erbe generator (Erbotom T) LiNA HF generator Adaptor – reusable – For Lamidey generator LiNA HF generator Adaptor – reusable – For Martin generator LiNA HF generator Adaptor – reusable – For Neomed generator LiNA HF generator Adaptor – reusable – For Siemens type 1 generator LiNA HF generator Adaptor – reusable – For Siemens type 2 generator	1 pc per box.
REC-617 (for 308)	Monopolar Silicone Cable - w. 4mm connector, 3m, reusable.	1 pc per box.

LiNA distributor in your country:



MD87093