

SIT 4

SIT 4 PLUS

SIT 4 PLUS E

The ENT
Examination and
Treatment Chair

An advantage for
doctor and patient



otopront®

SIT 4 SIT 4 PLUS SIT 4 PLUS E



ENT Examination and Treatment Chair

Use of high-quality materials and clear lines make this examination and treatment chair unique in both form and function.

Easy access to all adjusting elements as well as a close positioning to the patient is thus ensured.

Electric adjustment of the seat height by either 20 or 30 cm allows the doctor to work in an ergonomic position whilst sitting or standing. Using the version with 30 cm lifting range no additional seat accessories are needed for examining children.

While the backrest of the SIT 4 PLUS is adjusted, the armrests and leg rest simultaneously follow the chair's sequence of movement until it reaches a complete horizontal position thus ensuring perfect support and patient comfort. In case of a collapsing patient the backrest can quickly be brought into the -10° Trendelenburg position.

The patient chair is also available with electric backrest adjustment (SIT 4 PLUS E) and electric rotation (SIT 4 PLUS ER).

Details



A special feature of the SIT 4 is the **forward slanting telescopic column**, as a result of which the patient is automatically brought

towards the doctor when the chair is raised. This ensures even better access to the patient.


The version with **30 cm lifting range** ensures to use the chair for adults as well as for children without additional seat accessories.



The optional version with a **rotatable seat** with 90° locating mechanism facilitates rapid changes of position. The padded armrests are easily folded back out of the way to permit heavy-weight or wheelchair-bound patients access.



Version with **fixed footrest** (SIT 4).

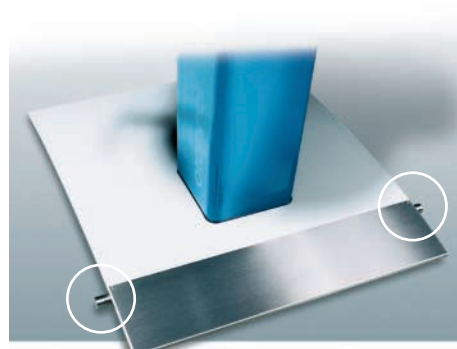

The height of the treatment chair is easily adjusted using the footswitches located on both sides of the base plate. There is also the option of an external footswitch. The wide stainless steel kick plate protects areas subject to particular wear.

SIT 4


SIT 4 PLUS

SIT 4 PLUS E

Technical Data

Lifting range	20 cm, optional 30 cm
Lifting mechanism	motor support
Seat height	51–71cm, optional 51–81cm
Rotation	360° with 90° locating mechanism, optionally electric rotation $\pm 180^\circ$
Locking brake	operated from the left and right side of the seat
Adjustment of the backrest	-10° forward to 10° beyond the horizontal position, markings at the 30° position, optionally with electric adjustment (SIT 4 PLUS E)
Armrests	simultaneous movement of armrests and leg rest synchronized by the movement of the backrest (SIT 4 PLUS)
Headrest	anatomically shaped, 15 cm adjustable height, completely removable
Casing	stable steel casing with stainless steel protection plate
Weight	100 kg
Power supply	110–230 V, 50/60 Hz, 550 VA
CE mark	accord. to MDD 93/42/EEC

All designs and specifications subject to change without notice.

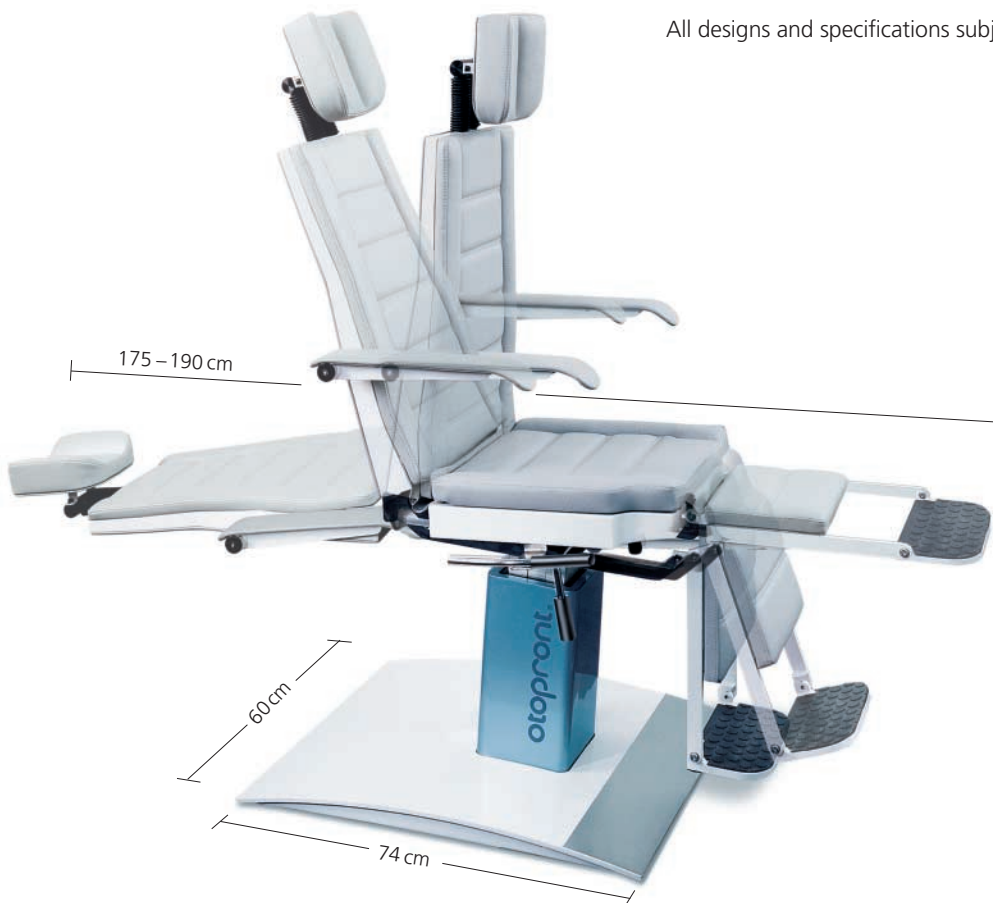

Happersberger otopront GmbH

Medical Equipment Manufacturer
D-65329 Hohenstein · P. O. Box · Germany
Tel. +49-6120-9217-0 · Fax +49-6120-921760
www.otopront.com · info@otopront.com

otopront®